
Quick Installation Guide

Follow the three steps below to unbox, assemble

① › ② › ③ ›

Digital Healthy Wheelchair Scale

（*Thank you for your choice*）

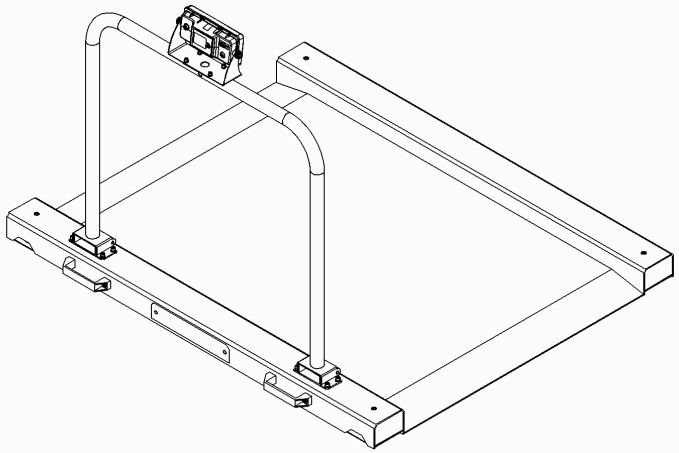

About LOCOSC PRECISION

LOCOSC PRECISION maintains a strong commitment to the ongoing development of precision weighing technology and products.

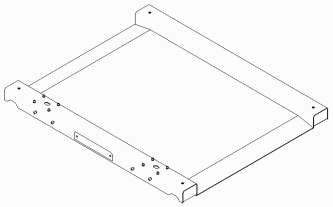
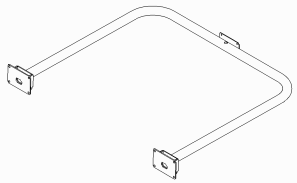
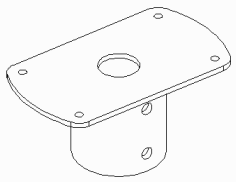
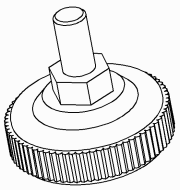
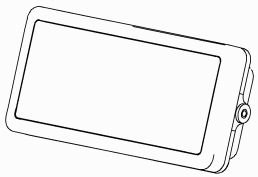
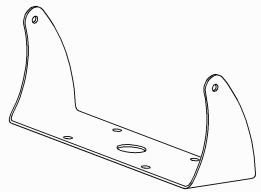
As a specialized manufacturer of precision weighing equipment, LOCOSC PRECISION serves customers across industrial, trade, and healthcare applications with a broad range of scales engineered for accuracy and reliability. This Wheelchair Weight Scale reflects that same standard of precision and quality control, so you can count on consistent, dependable measurements every time.

Table of Contents

1. Parts List
2. Accessory List
3. Assembly Instructions

1.Parts List






Before beginning assembly, unpack all components and check them against the list below.

A Base Platform	B Handrail	C Mounting Bracket (×1)	D Leveling Feet (×4)	E Weight Indicator (×1)	F U-Shaped Bracket (×1)
					

↓ CONTINUE TO STEP

2. Accessory List

The following hardware and accessories are included for assembly. Confirm the quantity of each item before starting.

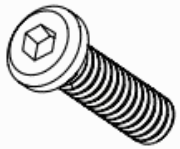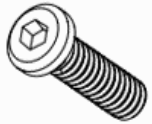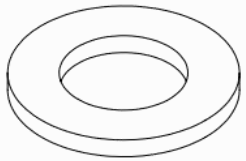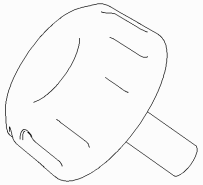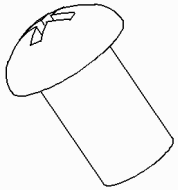1 Hex Socket Screw M8×16 (×4)	2 Hex Socket Screw M5×12 (×8)	3 Flat Washer (×8)	4 Spring Washer (×4)	5 5 mm Allen Wrench	6 6 mm Allen Wrench	7 Knob Screw (×2)	8 Phillips Pan Head Screw M5×8 (×4)
							

3.Assembly Instructions

Step 1 Insert the four feet into the corresponding holes at each corner, as shown.	**Step 2** Attach the pair of handles using M8×16 hex bolts.	**Step 3** Attach the handrail using M6×12 hex socket bolts.	**Step 4** Attach the U-shaped bracket to the handrail with the screws.	**Step 5** Mount the indicator to the bracket with screws, then connect the cable.
